

# Spex Standard Contour Cushion



## CLINICAL


Designed for the user who has moderate to severe postural needs and pressure requirements. Excellent choice for individuals who:

- Have changing postural needs over time
- Have a symmetrical or asymmetrical sitting posture that requires correction or control
- Require pelvic support and pressure relief/distribution
- Require a high degree of pelvic stability and specific contouring
- Need moderate control of abduction or adduction propensities

## FUNCTIONAL

- Created as the solution-solver for asymmetrical postures
- Strategic contouring pads offer complete customisations on-the-spot for a perfect fit for the client
- Gel-foam ischial well evenly distributes pressure under the pelvis
- Coccyx channel for enhanced pressure relief (for sizes 16" +)
- Dual-layer cover reduces shear forces, and breathable air-mesh on cushion sides counters moisture
- Inner incontinence cover included
- Customisable for many requirements such as leg length discrepancies

## CUSHION SIZE CHART

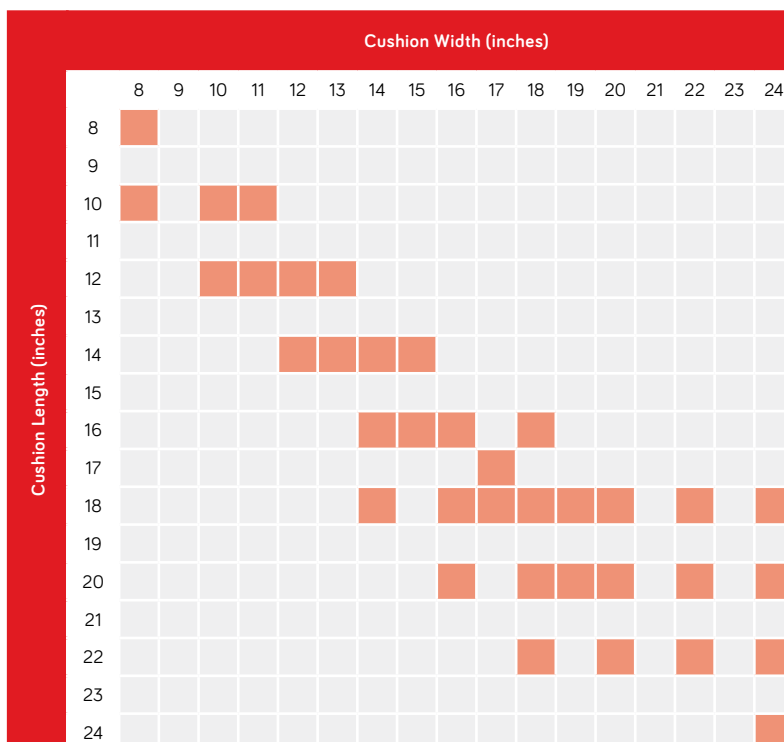

- Red squares are standard sizes.
- Non-standard or custom sizes are available on request.

## CUSHION DETAILS

Width	Height	Capacity	Trough Height	Product Weight
10"	2"	200lbs (90kg)	30mm	(16x16)
11"-15"	2.5"			<b>1.5kg</b>
16"-17"	2.75"	350lbs (160kg)		
18"-20"	3"			

## PRODUCT PERFORMANCE

Pressure Care	Red	Red	Red	4	
Positioning	Red	Red	Red	4	
Lateral Stability	Red	Red	3		
Offloading	Red	Red	Red	4	
Changing Needs	Red	Red	Red	4	

- Note: This is a guide only and ratings are dependent on correct use of product.

# Spex High Contour Cushion

The Spex High Contour Cushion offers enhanced thigh and abductor support. Key solution to positioning and levelling the pelvis correctly and effectively. Accomplish on-the-spot modifications with the cushion's strategic positioning base.

## CLINICAL

Designed for the user who has moderate to severe postural needs and pressure requirements. Excellent choice for individuals who:

- Need high control of abduction or adduction propensities
- Have changing postural needs over time
- Have a symmetrical or asymmetrical sitting posture that requires correction or control
- Require pelvic comfort and pressure relief/distribution
- Require a high degree of pelvic stability and specific contouring

## FUNCTIONAL

- Positioning base promotes anterior pelvic tilt and active reach to engage in activities
- Created as the solution-solver for asymmetrical postures
- Strategic contouring pads offer complete customisations on-the-spot for a perfect fit for the client
- Gel-foam ischial well evenly distributes pressure under the pelvis
- Coccyx channel for enhanced pressure relief (for sizes 16" +)
- Dual-layer cover reduces shear forces, and breathable air-mesh on cushion sides counters moisture
- Inner incontinence cover included
- Customisable for many requirements such as leg length discrepancies



Bespoke modifications can be scripted according to needs of the user.



## SPEX STRATEGIC POSITIONING BASE




Configure to suit many different postures of persons with disabilities or sitting issues. By adjusting the strategic positioning base, the cushion surface is contoured to improve overall posture, pelvic alignment, sitting ability, comfort and pressure reduction.




## PROMOTES ANTERIOR PELVIC TILT AND ACTIVE REACH TO ENGAGE IN ACTIVITIES

Cushion can be wedged to slope to front with cells added at rear and removed at front.

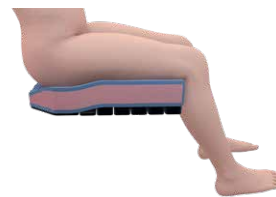

## MANAGEMENT & POSITIONING

Build cells up under thighs for enhanced pre-ischial support and control.

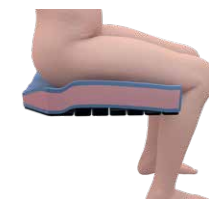

## PELVIC ROTATION

To help control pelvis from rotating, create a build-up under the affected leg (in front of the ischial).




Unique dual-layer **SHEAR REDUCTION** fabric surface

Accentuated thigh/adductor and abductor **CONTROL AND SUPPORT**

**CUSTOMISE** for all sitting postures

Maintain **PELVIC STABILITY** with adjustable positioning base

### CUSHION SIZE CHART

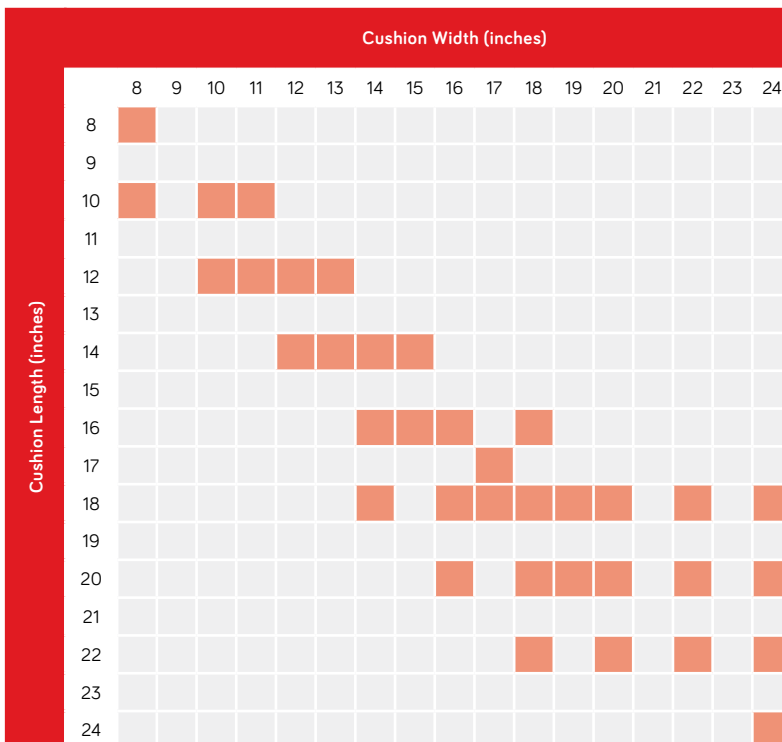

- Red squares are standard sizes.
- Non-standard or custom sizes are available on request.

### CUSHION DETAILS

Width	Height	Capacity	Trough Height	Product Weight
10"	2"	200lbs (90kg)	50mm	(16x16)
11"-15"	2.5"			
16"-17"	2.75"	350lbs (160kg)		<b>1.6kg</b>
18"-20"	3"			

### PRODUCT PERFORMANCE


Pressure Care	Red	Red	Red	4
Positioning	Red	Red	Red	4
Lateral Stability	Red	Red	Red	4
Offloading	Red	Red	Red	4
Changing Needs	Red	Red	Red	4

- Note: This is a guide only and ratings are dependent on correct use of product.

# Codes for Standard Sized Cushions

## Our cushion codes are meaningful!

Here's how to generate a product code to place your order.

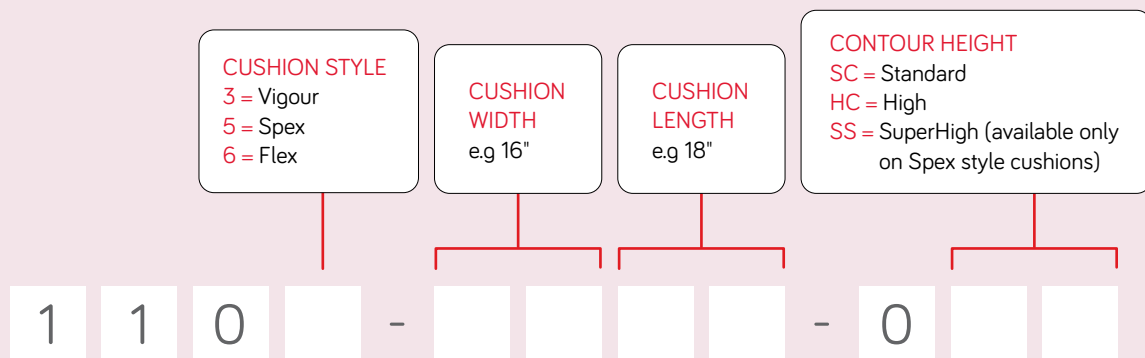


### EXAMPLE

A standard Contour Spex cushion of width 18" and length 20", would have the code 1105-1820-OSC.

## SPEX & VIGOUR CUSHION MEASUREMENTS GUIDE

This FRONT-VIEW of a Spex/Vigour cushion, shows the different measurements referred to on the subsequent product pages.

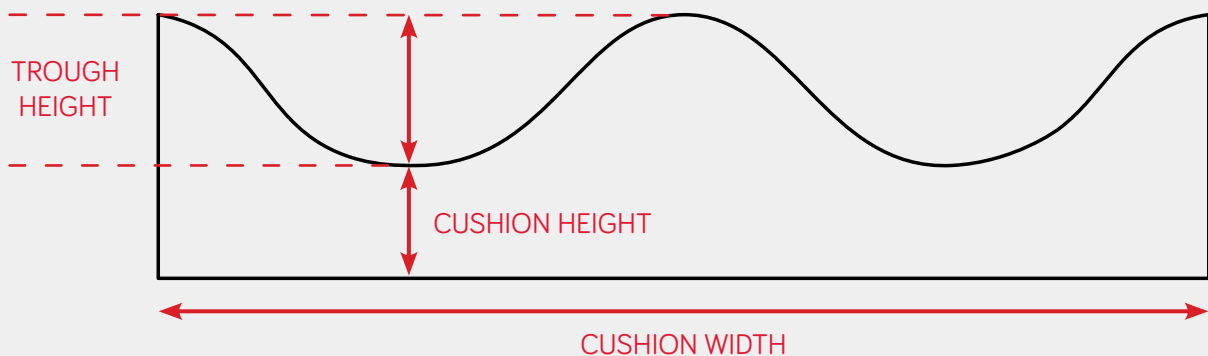

## Stand out from the crowd!

Upgrade your cushion with a splash of colour, for no extra cost!

We offer 8 colour options for our Spex cushion outer covers alternative to the standard black fabric. Simply add the upgrade codes to your cushion order.



**Granite Grey**  
1111-0491-090

**Ocean Blue**  
1111-0491-060

**Royal Purple**  
1111-0491-077

**Calypso Pink**  
1111-0491-028



**Chilli Red**  
1111-0491-022

**Mango Orange**  
1111-0491-030

**Bumblebee Yellow**  
1111-0491-044

**Leaf Green**  
1111-0491-054

## Spare cushion covers

All our cushions come standard with an inner incontinence cover and an outer cover, however if you would like a spare incontinence or outer cover, these are available.

See our scripting and ordering guide for ordering codes or get in touch with your local Spex distributor.

